
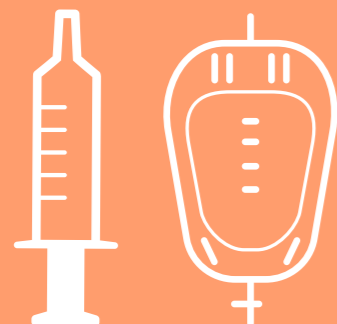


constantly creating better

ugofix.™
sleeve

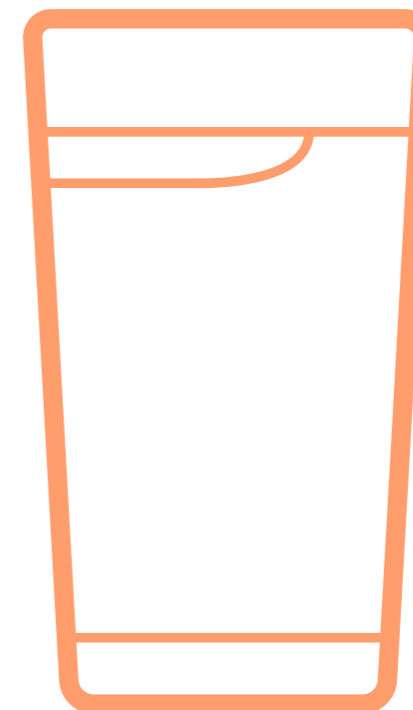

We're specialists in
urology products.

Find out more about our Ugo Urology and
OptiLube catheterisation jelly products by visiting
www.optimummedical.co.uk.



Optimum Medical Solutions Ltd
Tennant Hall, Blenheim Grove,
Leeds, LS2 9ET, United Kingdom
T: +44 (0) 845 643 5479
enquiries@optimummedical.co.uk
www.optimummedical.co.uk

EC REP

MT Promed Consulting GmbH
Altenhofstrasse 80
66386 St. Ingbert Germany



OMS-UFS-102-Apr20



ugofix sleeve

At **Optimum Medical** we understand the need for secure and reliable catheter and leg bag fixation.

We have developed the **Ugo Fix Sleeve** to provide comfortable support giving the patient the confidence to move around freely and safely.

The **Ugo Fix Sleeve** is designed to support leg bags attached to urethral catheters, suprapubic catheters and urinary sheath systems.

KEY STRENGTHS

- Washable, knitted, fibre blends provides strength and durability.
- Contours to leg shape for even distribution of urine and reduced noise.
- Graduated lines of concentrated elastane for support and controlled expansion when bag is filling.
- External seam to prevent skin damage and provide comfort.
- Strengthened outlet hole to maintain continuous support.
- Increased width and strength to top and bottom bands for support and fixation.
- No unnecessary tube holes in the top seam therefore preventing disconnection of the link drainage system.
- Colour coded top seam to assist with correct fitting and re ordering.

WARNING!

All patients should undergo an assessment for suitability and sizing prior to using a leg bag holder.

ORDERING INFORMATION

Before fitting, it is important that you measure around the widest part of your thigh or calf and use the correct size in the table below.

Colour	PR Code	PIP Code	Product Description	Inner Pack
Yellow	3005	394-0178	Ugo Fix Sleeve Small - 24cm-39cm	Box of 4
Blue	3006	394-0160	Ugo Fix Sleeve Medium - 36cm-55cm	Box of 4
Brown	3007	394-0152	Ugo Fix Sleeve Large - 40cm-70cm	Box of 4
Green	3008	394-0145	Ugo Fix Sleeve Extra Large - 65cm-90cm	Box of 4
	3027	409-6699	Ugo Fix Sleeve XXL - 75cm-105cm	Box of 4



We recommend that the **Ugo Fix Sleeve** is changed daily or sooner if the strap becomes soiled. To make this easy we have packaged the straps in sets of 5 to assist with the daily wear and wash cycle.

To maintain the **Ugo Fix Sleeve's** optimum performance, we recommend that the patient washes the strap in a laundry bag. To assist with this we provide a laundry bag in each box of **Ugo Fix Sleeves**.

THE IMPORTANCE OF CATHETER AND LEG BAG FIXATION

Proper securement is essential to prevent urethral trauma, pain and discomfort.

Wilson (2016) identifies 'To manage the risk of trauma and prolonged traction, adequate securement and fixation are required' and 'Securement refers to the use of leg-bag straps or mesh sleeves for supporting leg-bags while fixation indicates either the straps or adhesive devices that anchor the catheter.'

The **Ugo Fix Sleeve** is ideally complimented by our **Ugo Fix Catheter Strap** or **Ugo Fix Gentle** (catheter clip).

Find out more about our fixation range at www.optimummedical.co.uk

Reference:

Wilson M (2016) Urinary catheterisation in the community: Exploring challenges and solutions. *British Journal of Community Nursing* October 2016 Vol 21, No 10 pp 442-445